


DL•A electric operating table features with integrated multi-function can fulfill the needs of various operations.

Features

1. The bed surface adopts imported carbon plastic platethe, the table can slide which can fulfill the needs of C type arm whole body examination.
2. Panel control and remote control, with battery can fit for 50 sets operating.
3. The mattress is wrapped with imported artificial leather. High-density memory sponge is formed in one piece, with no seam.
4. Imported LINAK motor without noise and with long life-span.
5. The base seat and upright column covering are made of high quality 304 stainless steel which can be easily cleaned and sterilized.
6. Foot treadle brake device which is easy for moving operating table.
7. The leg part use spring gas control ,separate leg part

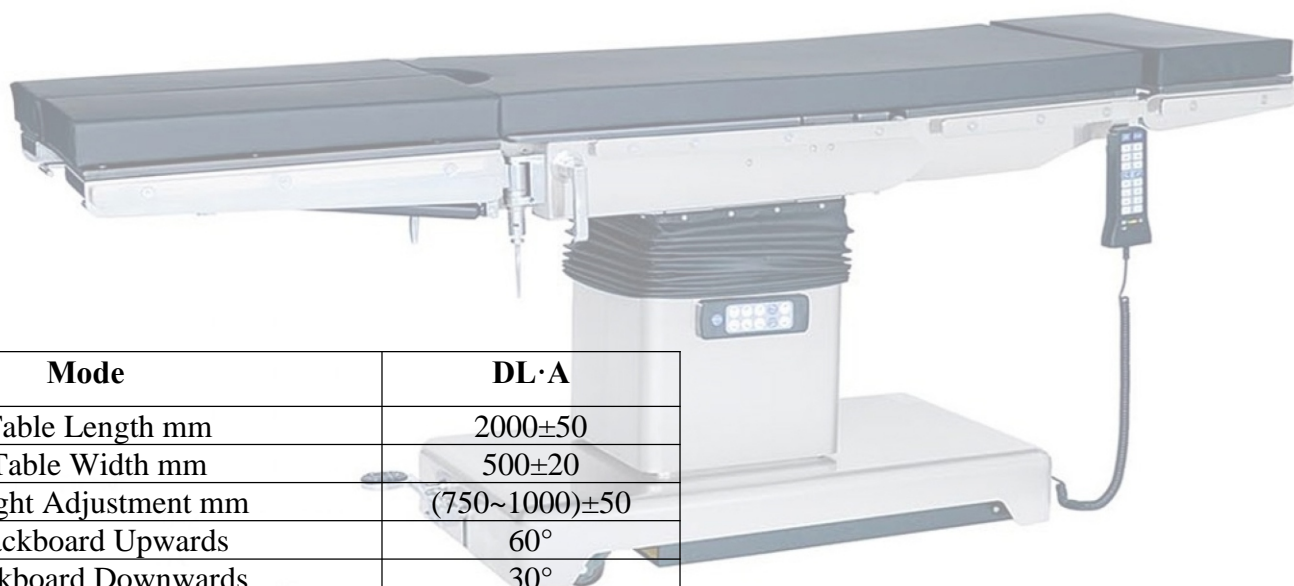

Standard Accessories: A pair of shoulder support frame, a pair of waist support frame, a pair of arm resting frame, a pair of leg support frame, a screen frame, a pair of pedal foot board, a waistband.

Mode	DL•A
Table Length mm	2000±50
Table Width mm	500±20
Height Adjustment mm	(750~1000)±50
Backboard Upwards	60°
Backboard Downwards	30°
Trendelenburg / Rev-Trend	20°
Tilt(Left) / Tilt(Right)	18°
Head board fold upward	40°
Head board fold downward	90°
Leg board fold upward / downward	20°/ 90°
Input power	300VA

Gross weight:306kg, Net Weight:263kg;
Volume:1300mm×750mm×1020mm.